



## SECTION 6.0: PARTS LIST

### 6.4 POWER SUPPLY ASSEMBLY

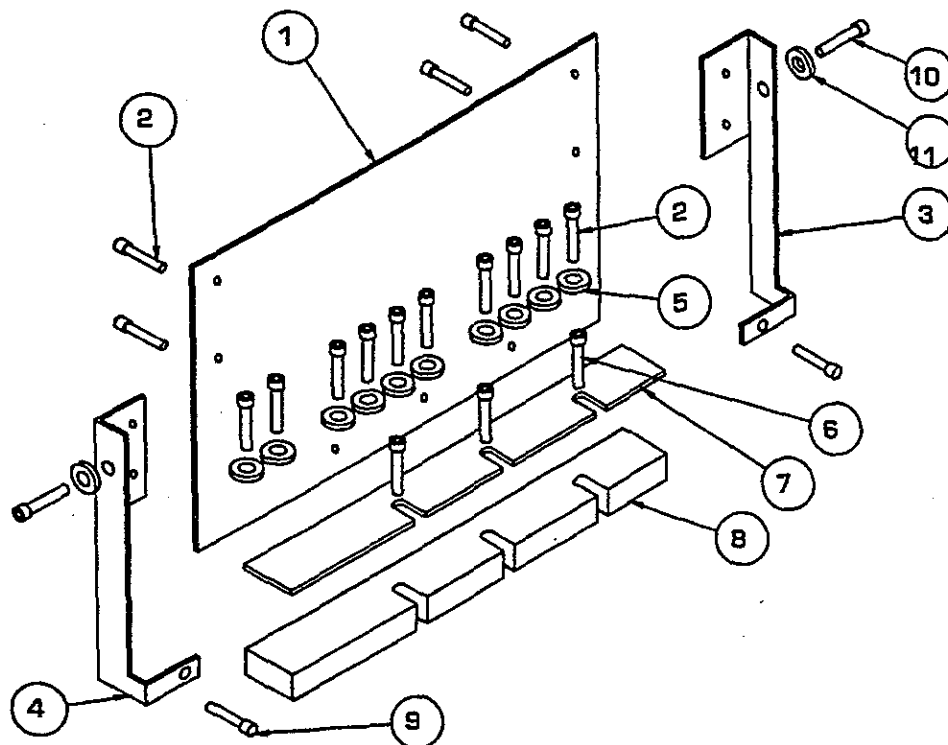

Item	P/N	Description	Qty
1	50160	Power Supply PC Board Assembly, V.I.P.	1
2	40013	Screw, 8 - 32X .500 HX CAP	14
3	20174R	Bracket, Mounting, Upright, Right	1
4	20174L	Bracket, Mounting, Upright, Left	1
5	45020	Washer, Insulator # 4 X .16	10
6	08833D	Screw, 8 - 32 X .375 SKTHDCP	3
7	20147	Insulator, Power Driver PC Board	1
8	20123	Bracket, Power Drive, PC Board	2
9	40018	Screw, 4 - 40 X .250 PNHD	2
10	03219D	Screw, 6-32 X .375 HXCAP	2
11	04383D	Washer, Nob INTL T LK WASH	2


## 6.5 EXHALATION VALVE

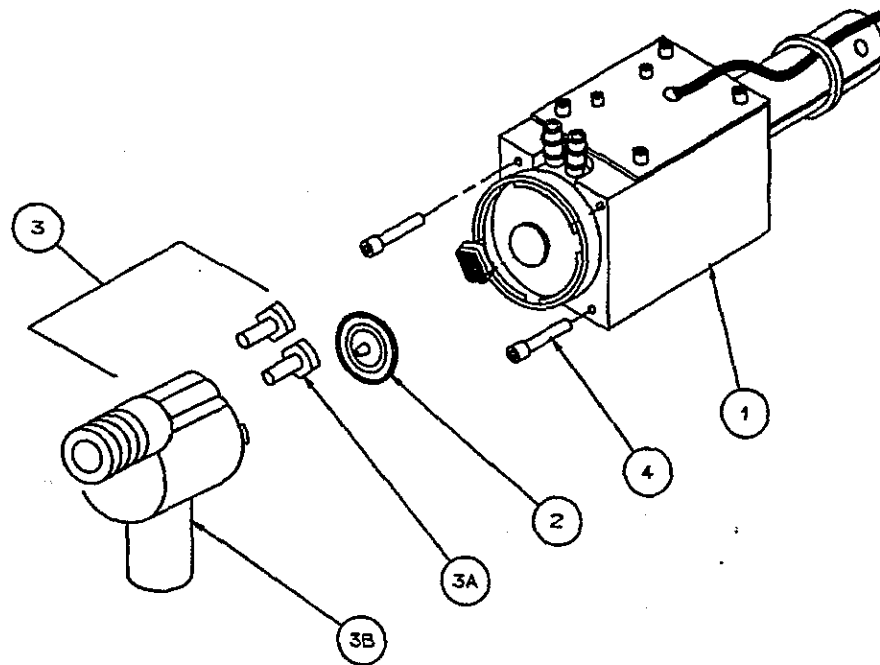

Item	P/N	Description	Qty
1	15029	Exhalation Valve Assembly	1
2	20107D	Diaphragm, L E	1
3	10088	Exhalation Valve Body Assembly	1
3A	20117	Seal, Barrel, L.E.V.	1
3B	10087	Exhalation Valve Body	1
4	03219	Screw, 6 - 32 X .375	2


## SECTION 6.0: PARTS LIST

### 6.6 MANIFOLD, V.I.P. BIRD ASSEMBLY

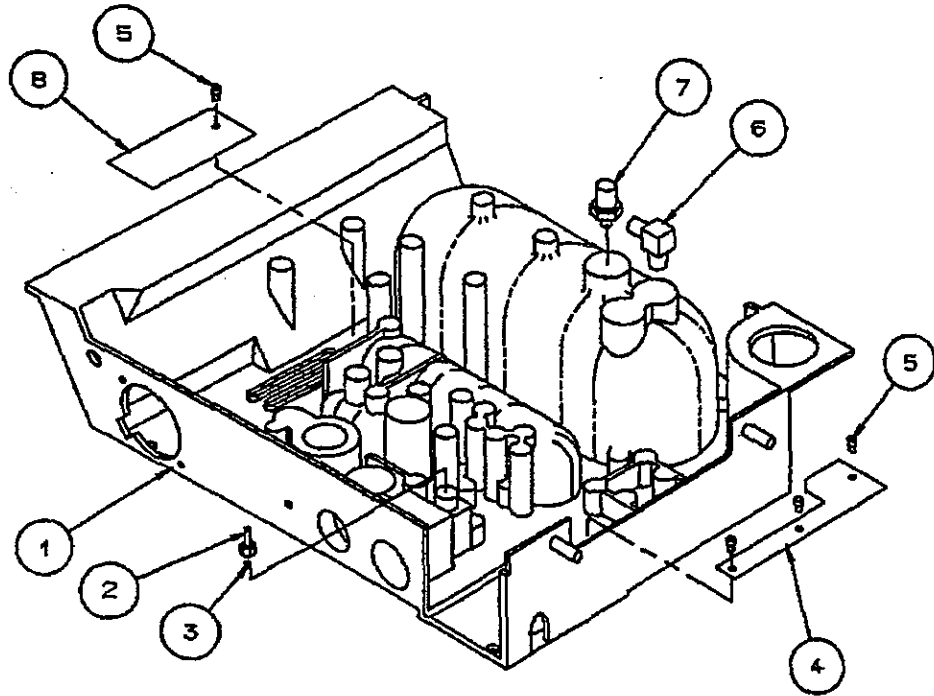

Item	P/N	Description	Qty
1	20125	Manifold, V.I.P Bird	1
2	00576D	Adapter, 10 x 32 x 1/8" Tubing	1
3	00114D	O-ring, .117 X .040	1
4	20293	Screen, Ventilator Right Hand Side	1
5	02865	Screw, 6-32 X .250 HXSKRDH	4
6	33615	Connector, 1/4 NPTM	1
7	09786	Valve, Pressure Relief	1
8	20292	Screen, Ventilator Left Hand Side	1


## 6.7 OUTLET PATIENT ASSEMBLY

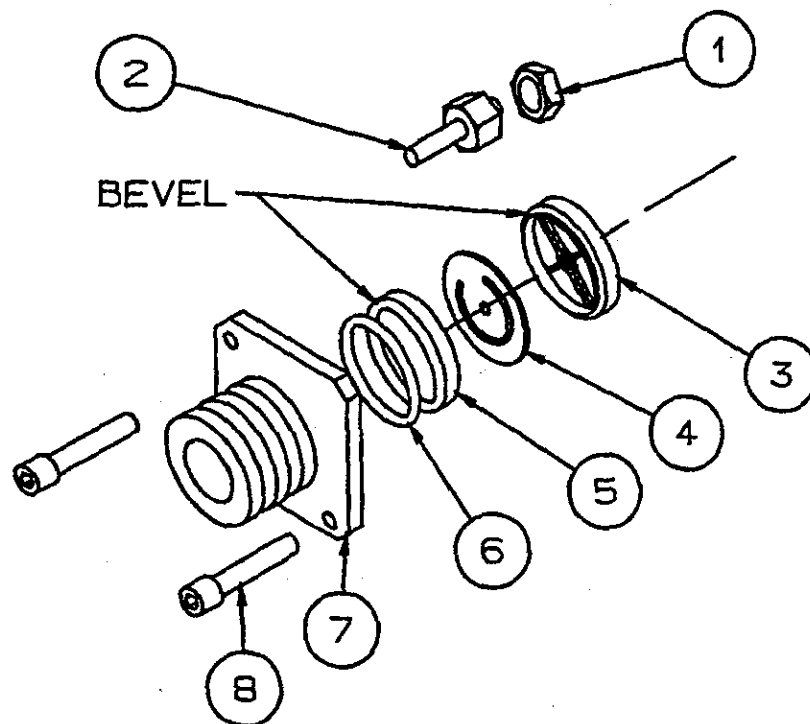

Item	P/N	Description	Qty
1	01537	Coupling, Panel 1/8 x 3/16 Tube	1
2	01056	Nut, 1/4 - 28 E-Keeps	1
3	05532	Retainer, Valve, Female	1
4	02168	Valve, Flapper	1
5	05531	Retainer, Valve, Male	1
6	02013	O-Ring, .864 x .070	1
7	10089	Outlet, Patient, V.I.P. BIRD	1
8	40013	Screw, 8 - 32 x 1.250 HXSCH	2

## SECTION 6.0: PARTS LIST

### 6.8 VALVE ASSEMBLY, RELIEF ADJUSTABLE

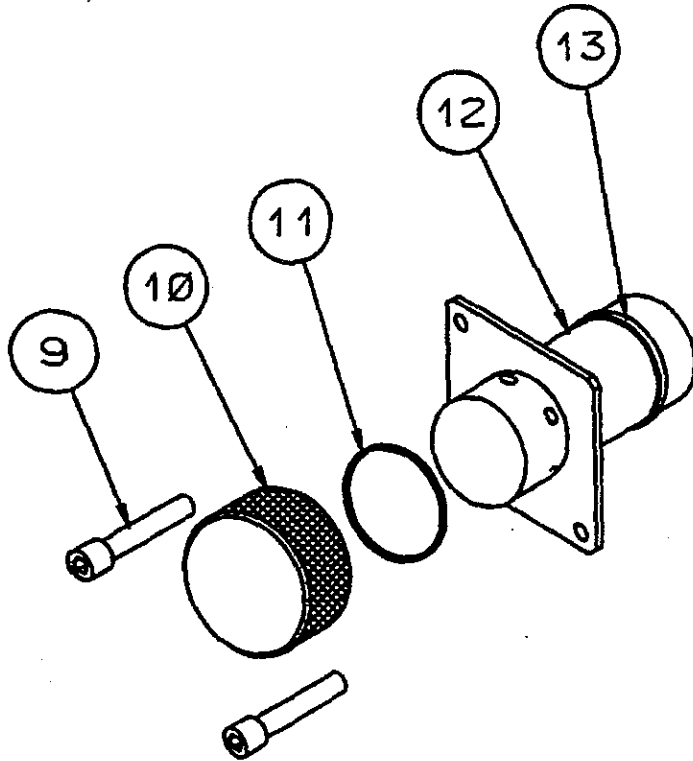


Item	P/N	Description	Qty
B	10092	Valve Assembly, Relief Adjustable	1
9	03219D	Screw, 6 - 32 X .375 HXCAP	2
10	10093	Knob Assembly, Relief Valve	1
11	30006	O-Ring, .989 X .070 Self Lube	2
12	20134	Body, Relief Valve	1

---

**THIS PAGE LEFT INTENTIONALLY BLANK**

### 6.9 SAFETY VALVE ASSEMBLY

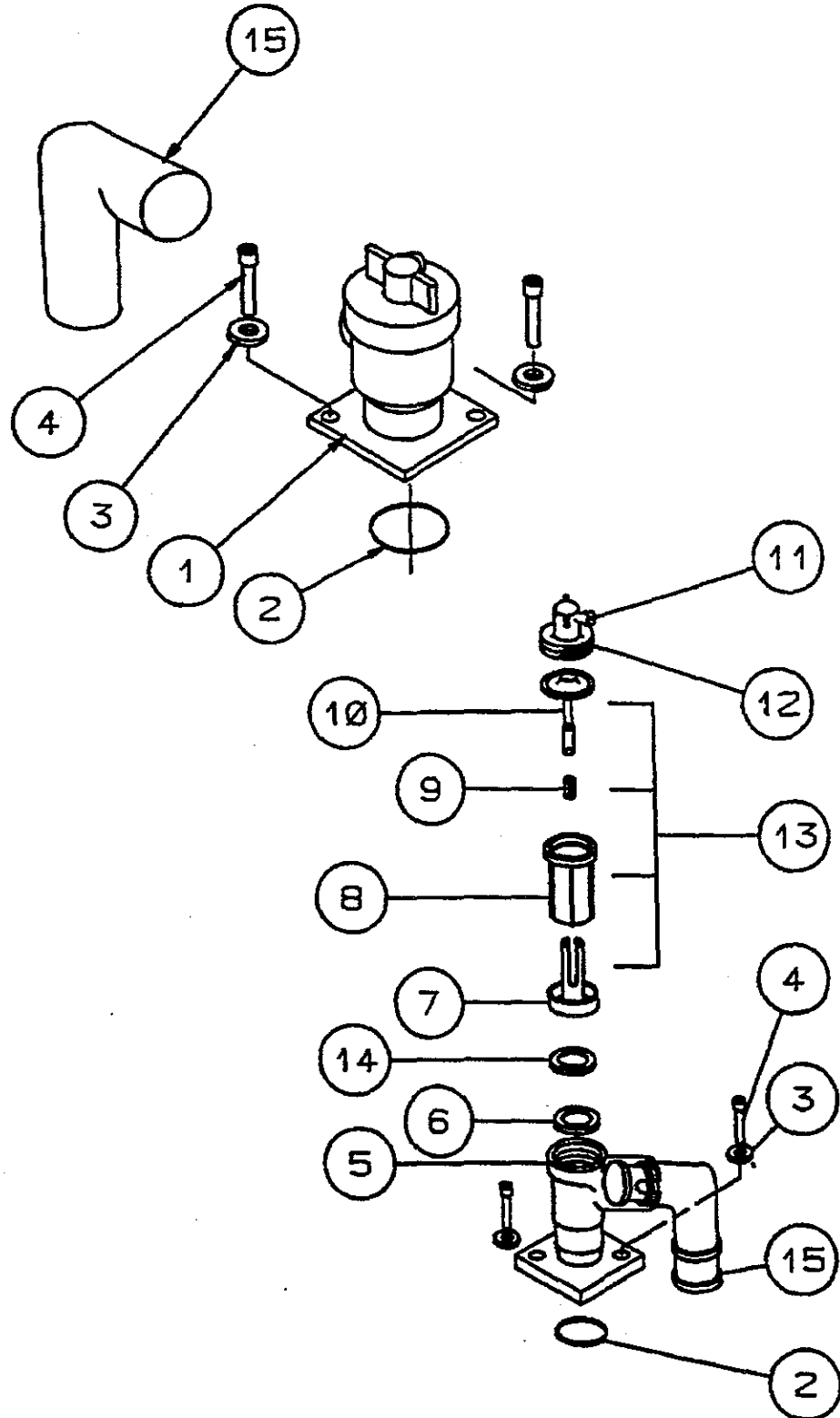

### 6.9 SAFETY VALVE - PARTS LIST

Item	P/N	Description	Qty
1	09754	Safety Valve Assembly	1
2	03373D	O-ring, .801 X 0.070 **	1
3	08964D	Washer, .148 X .272 X .030 **	2
4	03219D	Screw, 6-32 X .375 HXCAP **	2
5	10158	Valve, Bodyplate Assembly *	1
6	20244	Seal, Safety Valve *	1
7	02578D	Valve, Mushroom *	1
8	02577U	Insert, Exhalation Valve Cartridge *	1
9	20115	Spacer, Safety Valve *	1
10	08580D	Diaphragm, Exhalation Valve *	1
11	04143D	Retainer, Tube Connector	1
12	00172A	Cap, Exhalation Valve Assembly	1
13	10004	Cartridge, Safety Valve	1
14	04501	Washer, 1.100 X .900 X .20	1
15	09603	Connection, Elbow Safety Valve to	1
		Manifold, 6400ST **	


\*\* Not a part of Item 1

\* Subparts of Item 13, cartridge safety valve



## SECTION 6.0: PARTS LIST


### 6.10 PURGE VALVE ASSEMBLY

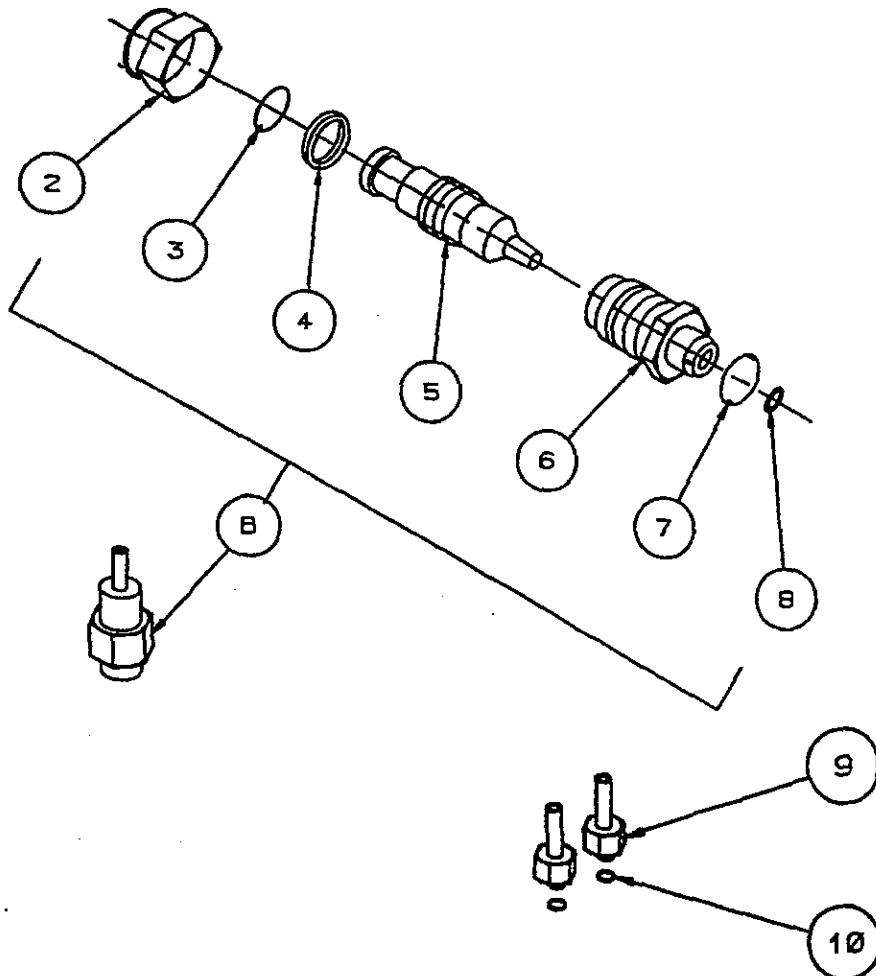

Item	P/N	Description	Qty
B	10180	Assembly, Purge Valve	1
2	00122	Nut, Packing Valve *	1
3	00109D	Washer, .193 X .281 X .031 *	1
4	00138D	O-Ring, .176 X .070 *	1
5	NA	Valve Stem, Flowrate *	
6	NA	Body, Valve*	
7	07849D	O-Ring, .313 X .051 *	1
8	00114D	O-Ring, .117 X .040 *	1
9	00576D	Adapter, 10-32/3A, 1/8" I.D. Tube	2
10	00114D	O-Ring, .117 X .040	2

\* Subparts of Item P/N 10180 Purge Valve Assembly

## 6.11 JET PUMP ASSEMBLY

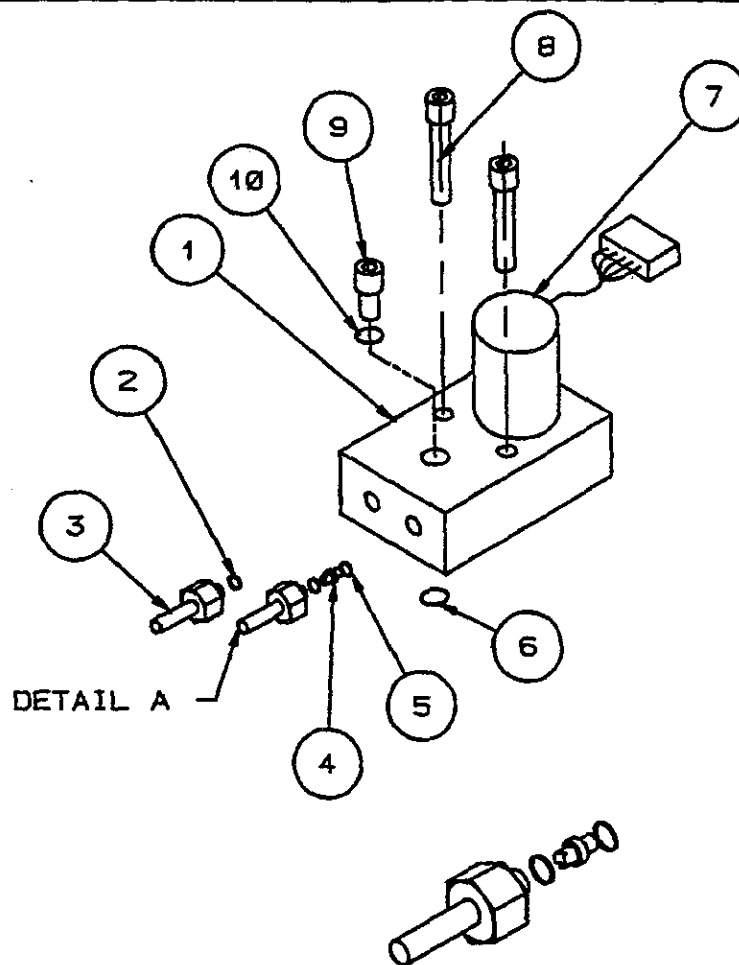


Item	P/N	Description	Qty
A	15009	Jet Pump Assembly	1
1	20142	Block, Jet Pump	
2	05999D	O-Ring, .180 X .040	3
3	04555	Adapter, 1/4 -28 X 1/2 Tube .104	2
4	08171	Orifice, .020 X 1/16 Tube Blue	1
5	00328D	O-Ring, .176 X .070	1
6	00138D	O-Ring, .070 X .070	1
7	15046	Solenoid Assembly, Jet Pump	
	0274D	O-Ring, .145 X .070	1
	03288	O-Ring, .551 X .070	1
8	40015	Screw, 8 - 32 Plug X 1.00 HXSCHD **	2
9	40026	Screw, 10 - 32 Plug, With Gasket	1
10	09538X	Tubing, 1/8 I.D., Flexible	6"

\*\* NOT PART OF (P/N 15009)

## SECTION 6.0: PARTS LIST


### 6.12 SAFETY VALVE SOLENOID ASSEMBLY

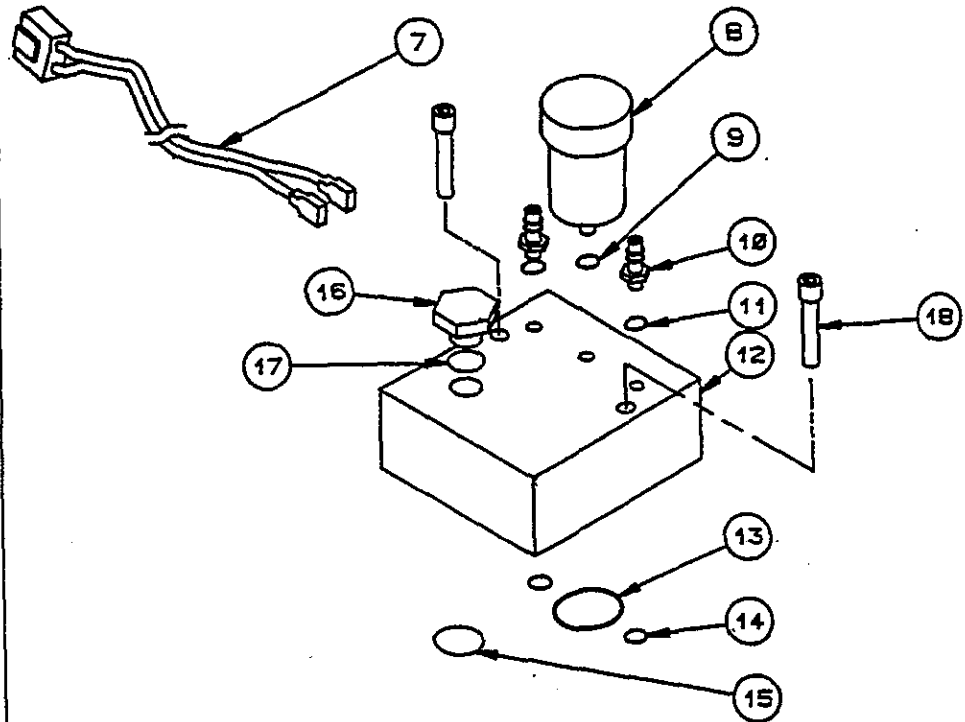

Item	P/N	Description	Qty
B	08940	Solenoid, Cable Assembly (W10)	
		Consists of	1
7	15120	Cable Assembly, Safety Solenoid Quick Disconnect	1
8	33603	Valve, High Flow, 3-Way	1
9	00274D	O-Ring, .117 X .070	1
10	00576	Adaptor, 10-32/3A, 1/8 Tube	2
11	00114D	O-Ring, .437 X .070	1
12	20145	Block, Solenoid Mounting	1
13	01943	O-Ring, .473 X .070	1
14	00138	O-Ring, .176 X .070	2
15	03375	O-Ring, .301 X .070	1
16	33527	Plug 9/16 - 18	1
17	03808D	O-Ring, .437 1.0 X .078	1
18	03217D	Screw, 8-32 X 1.375 HXCAP*	2

\* NOT PART OF (P/N 15013)


### 6.13 SIDE RAILS ASSEMBLY

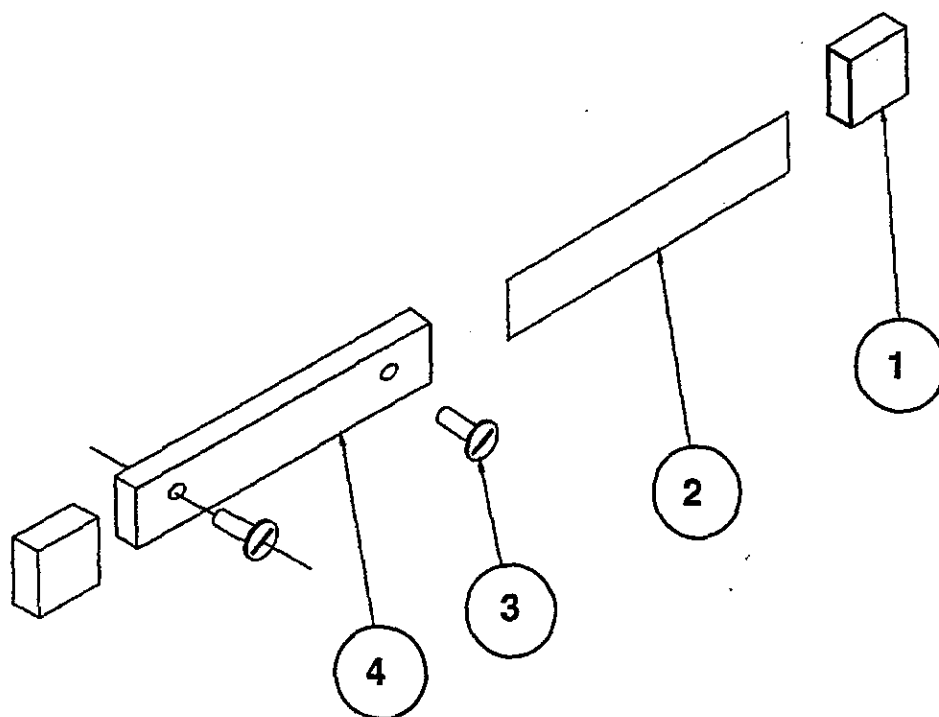

Item	P/N	Description	Qty
1	09342	CAP, End Rail	4
2	09391	Insert, Side Rail	2
3	03224	Screw, 1/4 - 20X .75 FLHSK	4
4	09696	Rail, Side	2


## SECTION 6.0: PARTS LIST

### 6.14 REGULATOR

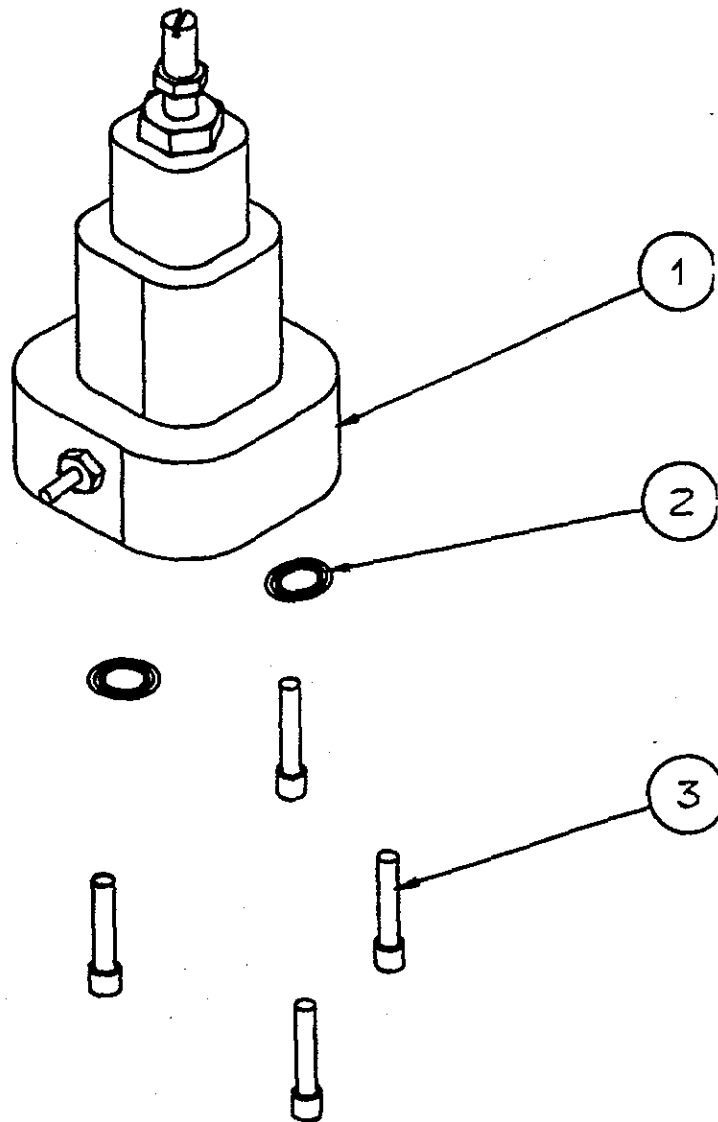

Item	P/N	Description	Qty
1	33566	Regulator, 2-25 PSI V.I.P.	1
2	03375D	O-Ring, .301 X .070 Silicone	2
3	03257	Screw, 8-32 X 1.5 HXCAP	4


## 6.15 CABLE ASSEMBLIES

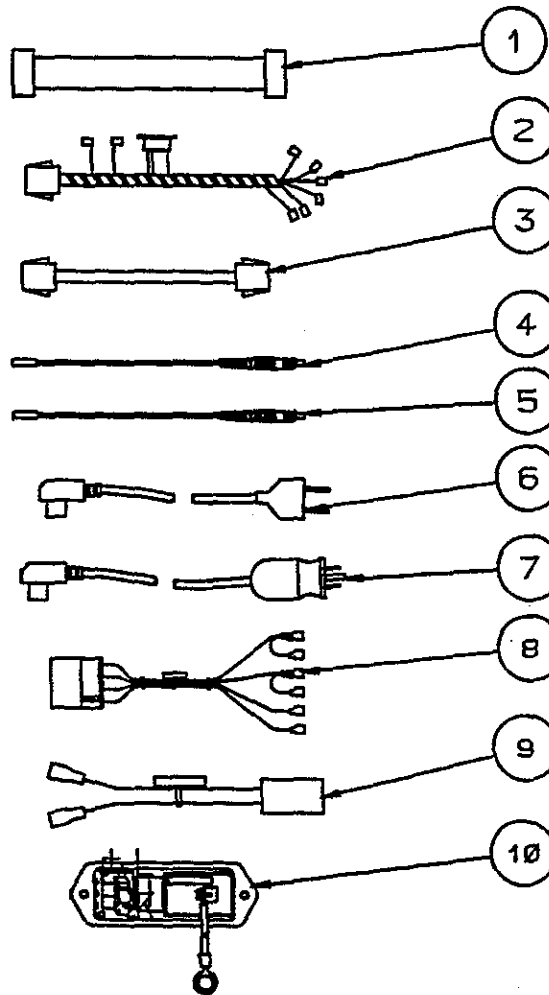

Item	P/N	Description	Qty
1	15044	Cable Assembly, Power Input V.I.P	1
2	15039	Cable Assembly, Power V.I.P.	1
3	15043	Cable Assembly, Data Interconnection	1
4	15016	Blue Cable, Fiber Optic 18 Inch	1
5	15126	Grey Cable, Fiber Optic 18 Inch	1
6	08925C	Powercord, 230V Right Angle	1
7	08925C	Powercord, 115V Right Angle	1
8	15116	Cable Assembly, Flow XDCR Manifold	1
9	15120	Cable Assembly, Safety Solenoid	1
10	15085	Cable Assembly, Power Entry W/O Switch	


## SECTION 6.0: PARTS LIST

### 6.16 MICROBLENDER

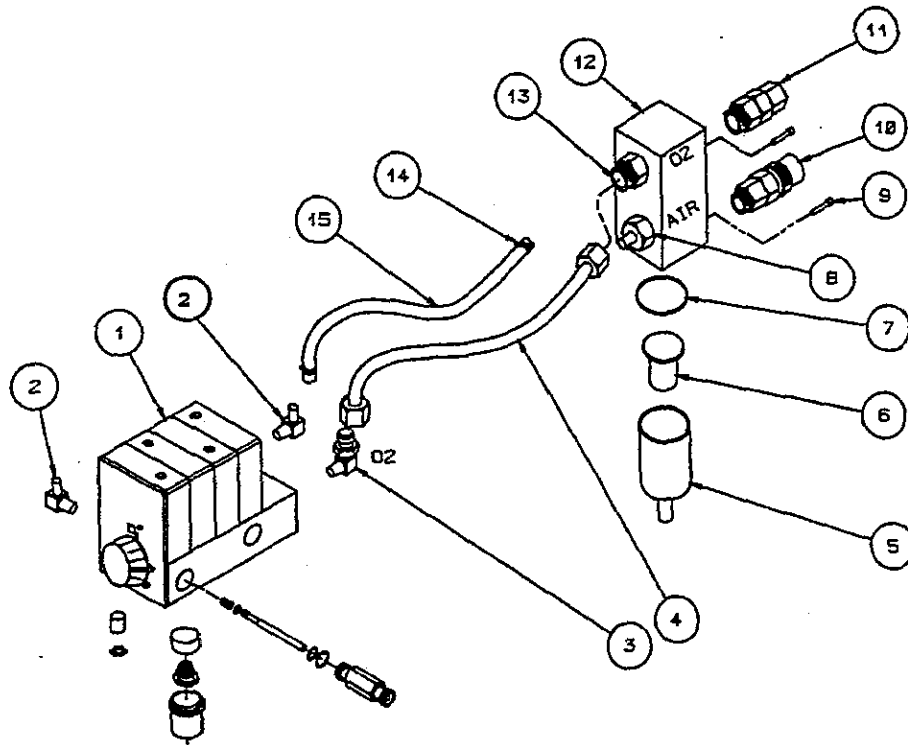

Item	P/N	Description	Qty
1	10160	MicroBlender	1
2	33521	Adapter, 90°, 1/8 NPT 1/4 Tube	2
3	00086	Adapter, 9/16 18 X 1/8 NPT ML	
		Connector, 1/8 NPT X 1/8 NPT MF EL	
4	09788X	Tubing 1/4 ID Reinforced	17'
5	09787	Clamp, Hose	
6	10197	Hose, Assembly, 0, .25 ID	
7	03864L	Inlet Assembly, 0, DISS	
8	10078	Inlet Assembly, Air	1
9	09835	Connector, 1/4 NPT - 1/4 Tube	
10	00610	ADAPTER, 9/16 18 X 1/4 NPT	
11	33519	Block, Air Inlet Filter	
12	06194	O-Ring, Filter	
13	06146	Filter, Element with Retainer	
14	06195	Housing, Filter	
15	05436	Alarm, Blender Bypass Assembly	
16	31000	Spring, .84 X .50 X .325 SS	
17	03903	Plug, Foam Alarm CAP	
18	03314	Ring, Star Retainer, Internal, .39	
19	03319	Muffler, Bleed	


## 6.17 FLOW CONTROL VALVE

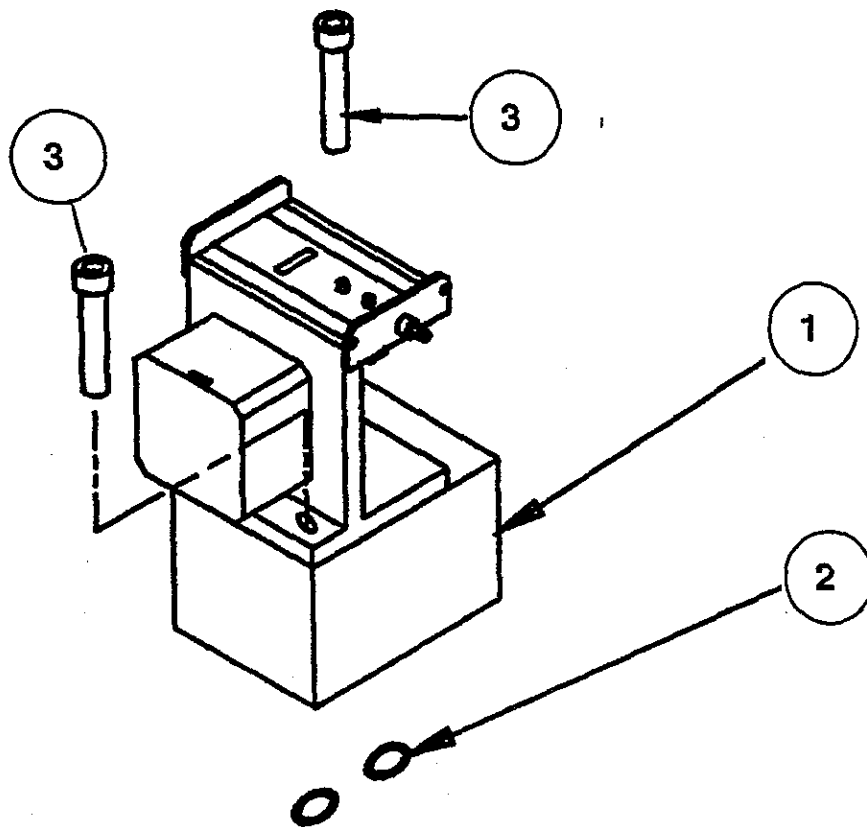

Item	P/N	Description	Qty
1	15005	Valve, HS/Flow V.I.P.	1
2	03374D	O-Ring, .364 X 0.070	2
3	03225	Screw, 8 - 32 X 2.00 HXCAP	2



### 6.18 TRANSFORMER ASSEMBLY

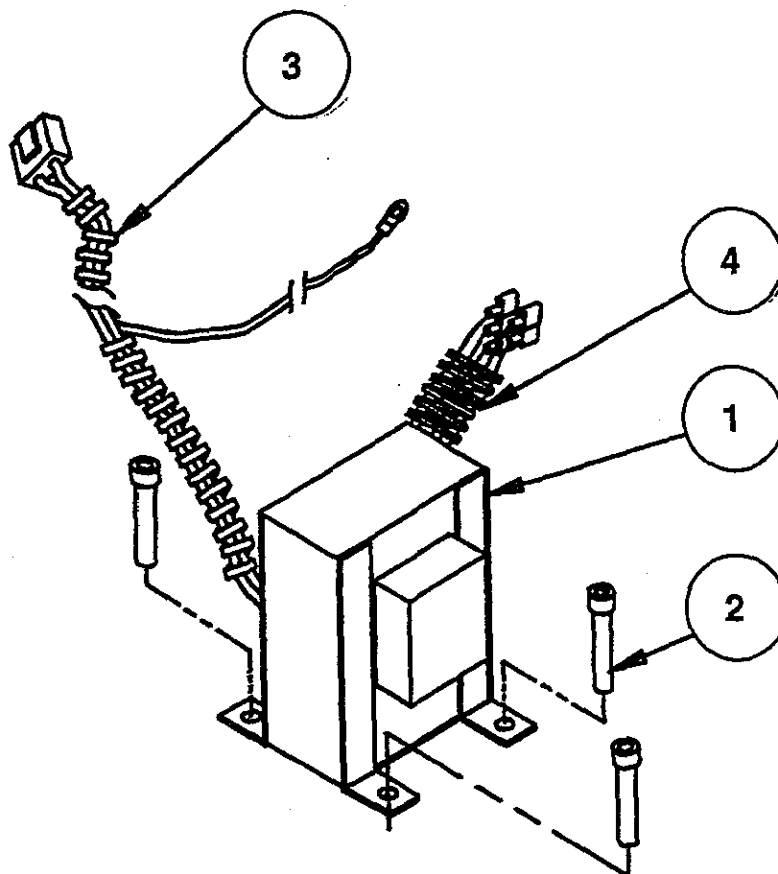


Item	P/N	Description	Qty
1	15114	Transformer with Cable Assembly (WG) IEC 60	1
2	08833D	Screw, 8 - 32 X .375 SKTHDCP	4
3	07803D	Tie Cable, Self Locking	1
4		Wrap, Spiral Nylon	12"

---

**THIS PAGE LEFT INTENTIONALLY BLANK**

## SECTION 6.0: PARTS LIST

### 6.19 PANEL ASSEMBLY, REAR


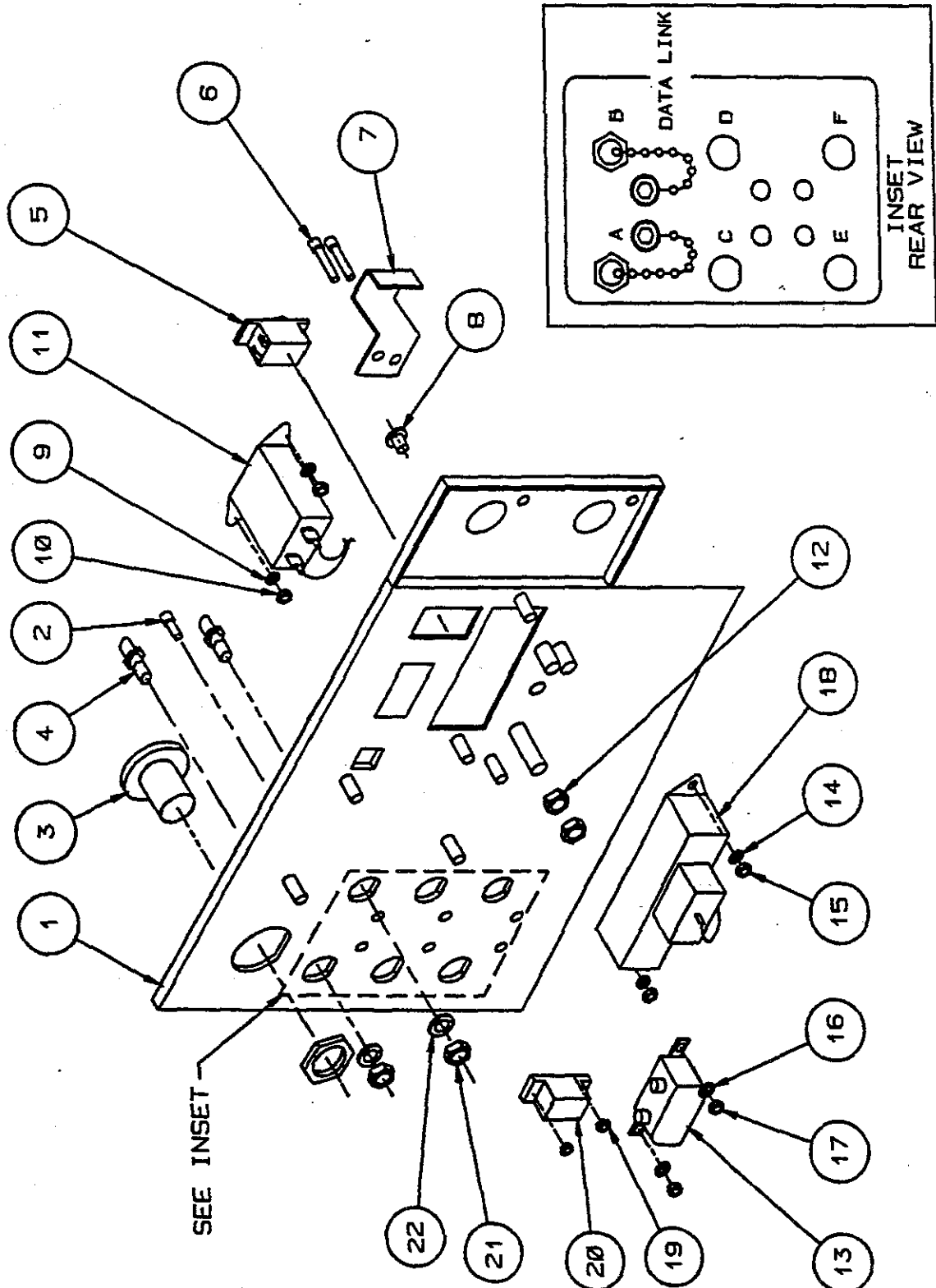


### 6.19 PANEL ASSEMBLY, REAR - PARTS LIST

Item	P/N	Description	Qty
1	20121	Panel, Rear V.I.P.	
2	03219	Screw, 6-32 X .375 HXCAP	6
3	03514	Connector Assembly, Male, 3 Pin (W1)	1
4	33628	Cap, Fiberoptic Connection	
5	68031	Switch, Voltage Regulator	
6	08833D	Screw, 8-32 X .375 SKTHDCP	2
7	06148	Bracket, Powercord	1
8	04647	Plug, Hole, .250 DIA Black	5
9	04383	Washer, No. 6 INTL T LK	12
10	07212	Nut, 6-32 X .093 Depth	6
11	09626	Meter, Hour 12 - 24 VDC	1
12	01066	Nut, 1/4 - 20 E KEPS	2
13	09625	Filter, AC Input	1
14	04383	Washer, No. 6 INTL T LK (Same as #9)	12
15	07212	Nut, 6-32 X .093 Depth (Same as #10)	6
16	07212	Nut, 6-32 X .093 Depth (Same as #10)	6
17	04383	Washer, No. 6 INTL T LK (Same as #9)	12
18	15085	Cable Assembly, Power Entry W/O SW	
19	03110	Nut, 4-40 UNC2B X .087 HEX	2
20	15044	Cable Assembly, Power Input	
21	33528	Connector, Bulkhead, Fiberoptic	2


## SECTION 6.0: PARTS LIST

### 6.20 TUBE ASSEMBLIES

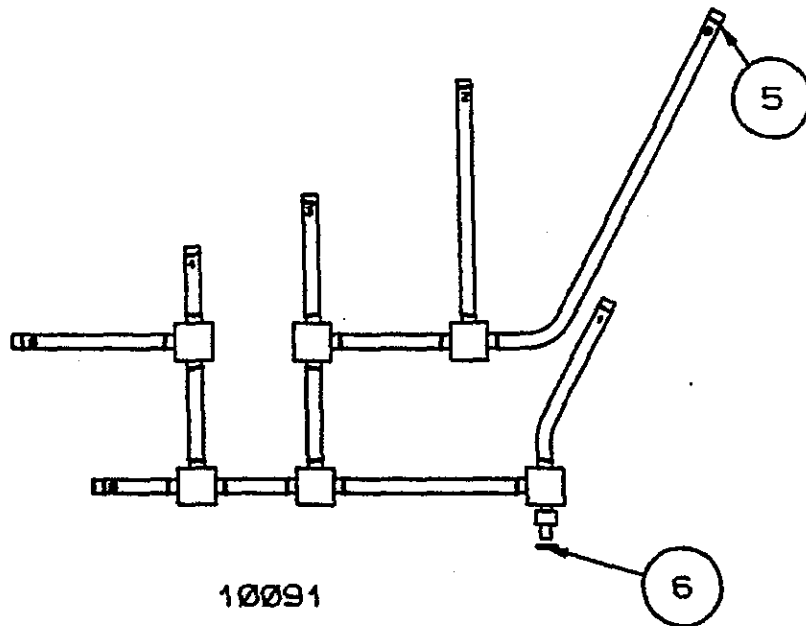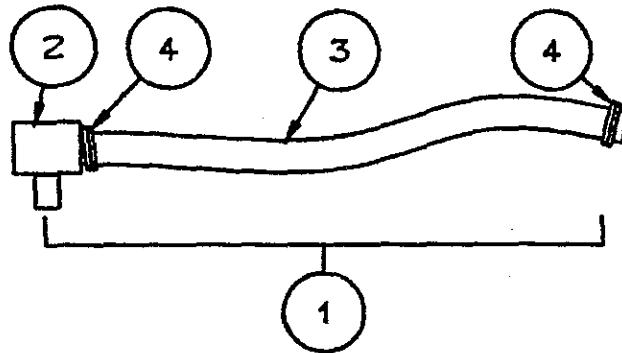


Item	P/N	Description	Qty
1	08934	Tube Assembly, Regulator Exhaust	
2	04506	Connector, 7.5mm X 1/8	1
3	09538X	Tubing, 1/8 I.D. Flexible	2
4	01587	Bracelet, 1/8 Tube	2
	10091	Tube Assembly, Prox	1
5	02189	Bracelet	24
6	06435	O-Ring, .097 X .026	1

## 6.20 TUBE ASSEMBLIES

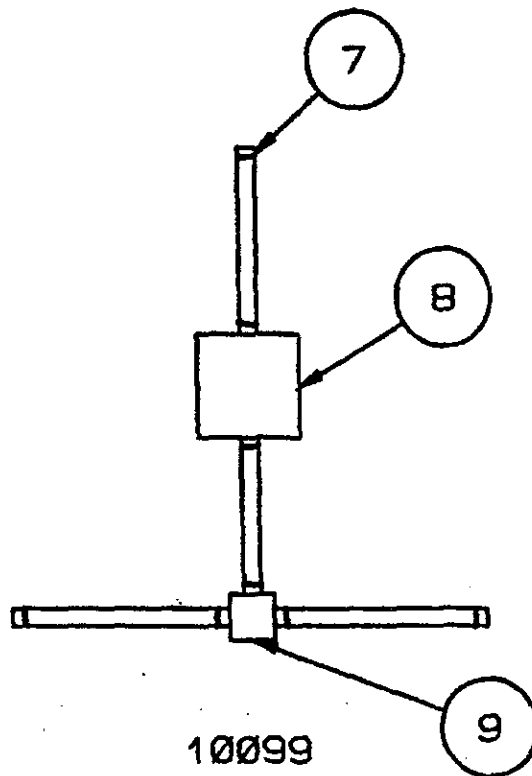

Item	P/N	Description	Qty
7	10099	Tubing Assembly, Muffler	1
7	02189	Bracelet, 5/16 Tube	1
8	15119	Muffler Assembly	1
9	00358	Connector, Tee 4.5mm X 1/8	1
10	08171	Orifice, Fixed .020 X 1/16 Tube Blue	1


6.21 EXPLODED VIEW, V.I.P. BIRD

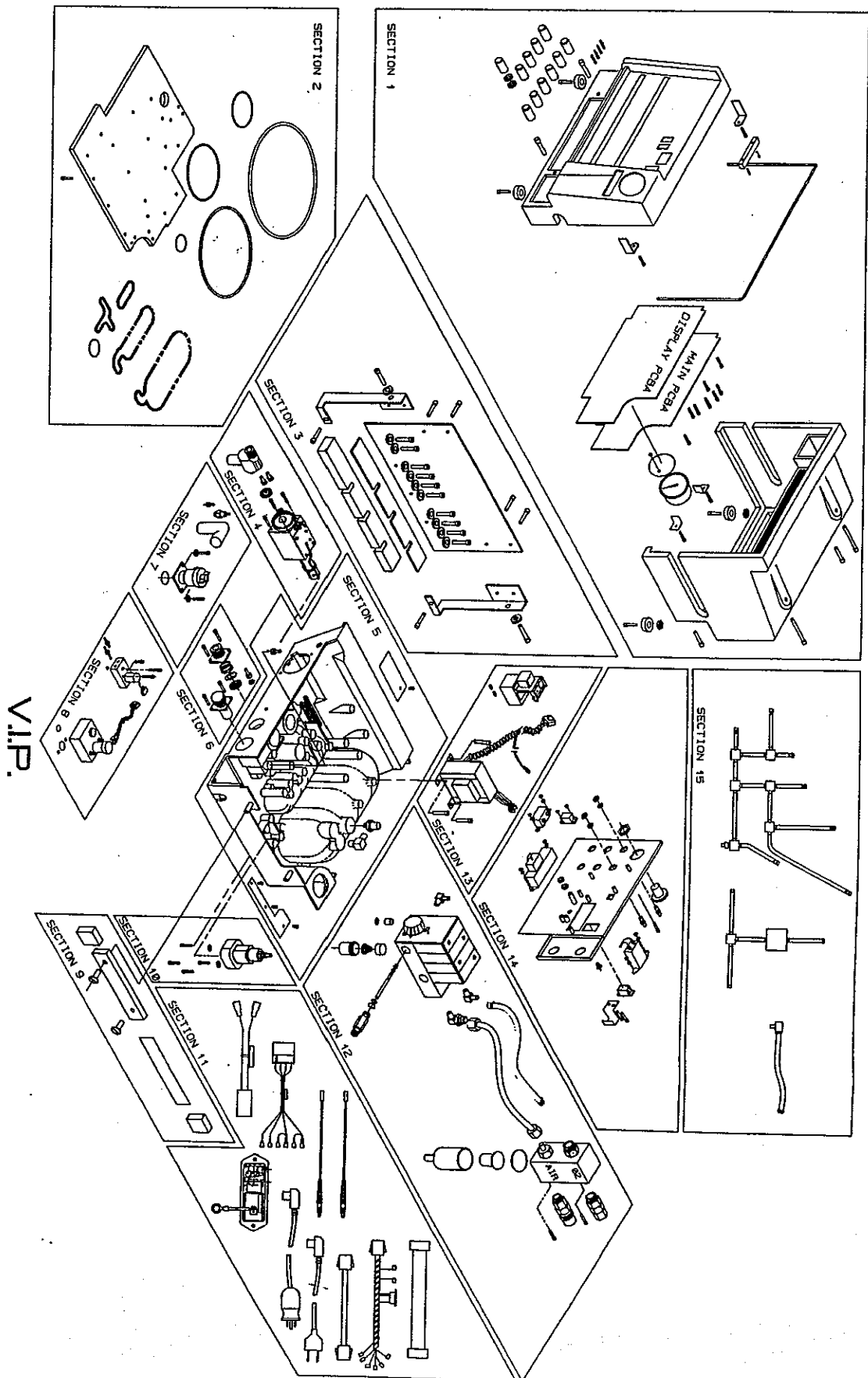

---

# Glossary

V.I.P. BIRD® Infant-Pediatric Ventilator  
Service Manual

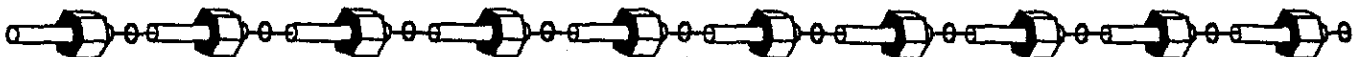



A/D .....	Analog to Digital
amp .....	Ampere
ANSI .....	American National Standards Institute
ASTM .....	American Society for Testing and Materials
BPM .....	Breaths Per Minute
C .....	Compliance
°C .....	Degrees Centigrade
CIRC .....	Circuit
cm .....	Centimeter
cmH <sub>2</sub> O .....	Centimeter of Water Pressure
cmH <sub>2</sub> O/LPS .....	Centimeter of Water Pressure Per Liter Per Second
CMV .....	Controlled Mechanical Ventilation
CPAP .....	Continuous Positive Airway Pressure
D/A .....	Digital to Analog
DISS .....	Diameter Index Safety System
DVM .....	Digital Volt Meter
ETO .....	Ethylene Oxide
°F .....	Degrees Fahrenheit
FIO <sub>2</sub> .....	Fractional Concentration of Inspired Oxygen
H .....	Height
Hz .....	Hertz
I.D. ....	Internal Diameter
I:E .....	Inspiratory Time to Expiratory Time Ratio
(I) .....	ON
IMV .....	Intermittent Mandatory Ventilation
kg .....	Kilogram
kg/cm <sup>2</sup> .....	Kilograms Per Square Centimeter
lb .....	Pound
LED .....	Light Emitting Diode
LPM .....	Liters Per Minute
ml .....	Milliliter
msec .....	Millisecond
mV .....	Millivolt
(O) .....	OFF
O.D. ....	Outside Diameter



## GLOSSARY OF ABBREVIATIONS


---

O <sub>2</sub> .....	Oxygen
PEEP .....	Positive End Expiratory Pressure
P/N .....	Part Number
PSI .....	Pounds Per Square Inch
PSIG .....	Pounds Per Square Inch Gauge
R .....	Resistance
SIMV .....	Synchronized Intermittent Mandatory Ventilation
t .....	Time
TCPL .....	Time Cycled, Pressure Limited
UL .....	Underwriters Laboratory
VAC .....	Volts Alternating Current
VDC .....	Volts Direct Current
Vt .....	Tidal Volume
W .....	Weight

---

# Bulletins

V.I.P. BIRD® Infant-Pediatric Ventilator  
Service Manual

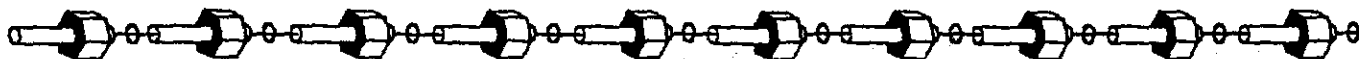

# 1994

# SERVICE BULLETIN

No: SE0003-94

TO: All VIP Service/Operator Manual Holders

SUBJECT: VIP Bird - Reported Vent Inops

FROM: Technical Service      DATE: April 22, 1994

## BACKGROUND

Bird Products Corporation is notifying users of the V.I.P-Bird Infant Ventilator that a cause for activation of the ventilator inoperative alarm could be one of the following:

1. Excessive (greater than twice the limit specified in the Food and Drug Administration's Medical Device Standard MDS-201-0004) electromagnetic interference.
2. Excessive Electrostatic Discharge between the device and operator.

The Operator's Manual and Service Manual both contain the following warning:

*The function of this unit may be adversely affected by the operation of equipment such as high frequency surgical (diathermy) equipment, defibrillators or short-wave therapy equipment in the vicinity.*

The end-user should take appropriate precautions to avoid close proximity to sources of Electrostatic Discharge or Radio Frequency Interference.

The V.I.P Bird meets the MDS-201-0004 Electromagnetic Compatibility Standard. The frequency of Vent Inop occurrences due to the effects of Electrostatic Discharge and Radio Frequency Interference is rare and appears to be isolated to only a few users. As such, Bird Products Corporation does not believe the population of installed ventilators to be at any significant risk in terms of safety or efficacy.

An enhancement kit (P/N 10326-Standard V.I.P.; P/N 10327-V.I.P. w/ Flow Sync) is available to users having ventilators that exhibit Vent Inop due to the influences from Electrostatic Discharge and or Radio Frequency Interference. This kit is already incorporated in all VIP-Bird infant ventilators with serial numbers above KDJ01354.

If you have any questions or concerns regarding this matter, please contact Bird Products Corporation Technical Service Department at (800) 934-2473.

